









Experience the Power of Precision in Your Daily Laboratory Routine



Hold Pipettes Firmly with a Universal Silicon Adapter

MICROLIT EASYFILL has a Universal Silicon Adapter, with internal knurling, that can easily accommodate glass/plastic pipettes from 1 ml to 100 ml and hold them firmly during liquid handling operations.



Perform Sensitive Laboratory

electronic pipette filling instrument,

Designed by our in-house team of product design engineers with

Single Knob Dual Dispensing Mode.

performing sensitive operations with

high precision, accuracy and reliability in practical laboratory environments.

ergonomics and intuitive handling in mind, it offers

MICROLIT EASYFILL helps in

is a simple yet powerful blend of sophisticated features

Operations with Assurance

MICROLIT EASYFILL, an

and functionality.

Forget Contamination with In-line Membrane Filter

A PTFE Hydrophobic In-line Membrane Filter (0.2 microns) prevents aerosol contamination and accidental liquid entry.



Know the Battery Status in One Go

Red LED flash denotes 'low battery' and green LED flash denotes 'charging'.



Switch Modes with a Single Knob

The instrument offers Dual Dispensing Mode with the help of a Single Knob.

For Gravity Mode : Press the knob halfway.

For Blow-out Mode : Press the knob fully.



Work with Ease

MICROLIT EASYFILL has a comfortable grip and exhibits high ease of use.



Adjust the Pump Speed During Liquid Handling

The instrument offers continuous pump speed adjustment during liquid handling operations.



Our Product Excellence





Our products reflect the design philosophy of our organization, combining the best performance parameters with the best user experience.

We were honoured by CII Industrial Innovation Awards and adjudged a place in the top 25 Most Innovative Companies of 2021.

We are also accredited by a number of global organizations and have earned worldwide credibility.

SPECIFICATIONS AND ORDERING INFORMATION











Switch Modes with a Single Knob

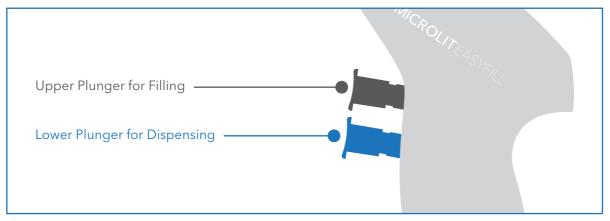


The instrument offers Dual Dispensing Mode with the help of a Single Knob.

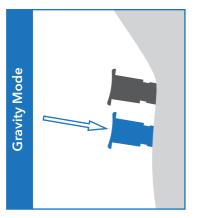




For Gravity Mode - Press the knob halfway.

This mode does not engage the motor and allows the user to dispense dropwise and slowly.

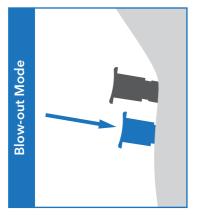

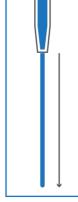

For Blow-out Mode - Press the knob fully. This mode engages the motor and

This mode engages the motor and dispenses the reagent with maximum speed.

Available in Four Colour Combinations

The instrument is available in four colour combinations.

*Refer to the specifications and ordering information.









MICROLIT

OUR REACH



YOUR SATISFACTION IS OUR PRIORITY

Microlit products follow GMP norms and are manufactured under comprehensive QA (Quality Assurance) and QC (Quality Control). However, if you are dissatisfied with the operation of any of our products, get in touch with us or call your nearest Microlit dealer for free replacement.

CONTACT US

Microlit INDIA

629 Pakramau, Kursi Road, Lucknow 226026, India Phone: +91 9918625629, Email: info@microlit.com www.microlit.com

Microlit LATAM

Business Place Ibirapuera, Alameda dos Jurupis,1005, Indianópolis, São Paulo, Brasil Phone: +55 11 93269-5932, Email: info-latam@microlit.com

Microlit USA

33 Wood Avenue South, Suite 600, Iselin, NJ 08830, USA Phone: +1 732 321 0852, Email: info-usa@microlit.com www.microlit.us